

Technical Sheet

Rev. 04 - 24/04/2015 - Mod. 72G2 - EN

POCKET 0,4 Lt.



GENERAL DATA

Disposable container with double opening for sharp and cutting objects disposal Item designed for low volumes of waste (care tray, trolley, domestic care, doctor's office)

- ▶ **Materiale:** Ecological plastic material, polypropylene copolymer coloured with masterbatch without heavy metals or cadmium Material compliant to REACH
- ▶ **Colour:** Yellow body and red lid in order to signal the dangerous content.
- ▶ **Thickness:**
 - Bottom: 2,26 mm
 - Side walls : 2,19 mm
 - Lid: 2,28 mm
 - Cap: 2,16 mm
- ▶ **Weight:** 105 gr.
- ▶ **Max load:** 0,26 Kg.
- ▶ **Nominal capacity:** 0,400 Lt.
- ▶ **Dimensioni esterne:** 124x63x165 mm H
- ▶ **Pallet dimensions:** 86x39 mm H



LOGISTIC DATA

- ▶ **Pcs. per carton:** 70 pcs. (body and lid not assembled)
- ▶ **Carton dimensions:** 60x40x36 cm H
- ▶ **Pcs per pallet:** 1680 pcs.
- ▶ **Pallet dimensions:** 120x80x228 H cm

CERTIFICATIONS:

- ▶ MARQUE NF
- ▶ Kite Mak - Compliance to ISO 23907
- ▶ Compliant to ADR - UN Certificate
- ▶ TRBA 250
- ▶ FDA K071517

INSTRUCTIONS

► Assembly:

The container is sold with body and lid assembled (ultrasounds assemblage)



► Temporary Closure:

The two caps on the lid can be closed in temporary position pressing with a finger.



At this stage, the cap can be reopened pressing its lower side

► Permanent Closure:


From the temporary closure position, firmly press on the "Push" symbol until it clicks.



At this stage, closure is permanent

► Needles removal:

Use the "V" removal device to separate needle and syringe

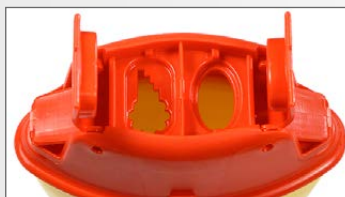

Use the "marguerite" device to unscrew the insulin pen



HOLDERS:

The container is provided with a handle



USE:

- Close the cap in temporary closure position
- Close permanently once full

ADVANTAGES:

- ▶ 2 openings
- ▶ Safe, hermetinc closure, easy to use and to handle
- ▶ Sharp and cutting object resistant
- ▶ Handle available to ease transport
- ▶ Insulin pens removal device
- ▶ Autoclavable at 134° (after this process, the container mantaines the features stated by NFX 30-500)
- ▶ Maximum filling level indicated
- ▶ Irreversible permanent closure
- ▶ To be disposed through incineration, releasing no toxic emissions

SAFETY INFORMATION:

- ▶ Do not use the product before fitting the lid
- ▶ Do not attempt to rescue objects already in the container
- ▶ Respect the filling level indicated: 80% of real capacity
- ▶ Use the container only with the appropriate lid
- ▶ Do not try to reopen the container once permanently and irreversibly sealed
- ▶ Follow the instructions of use to ensure the safety of the container
- ▶ Container to dispose according to the existing law
- ▶ Use handles to move the container