



"CICO"
can't intubate,
can't oxygenate


ScalpelCric

A surgical set for cricothyrotomy according to the scalpel technique


/ ScalpelCric

Use


ScalpelCric has been developed to match the latest recommendations from the DAS (Difficult Airway Society) 2015 guidelines. Based on the simple description of "**stab, twist, bougie, tube**" the ScalpelCric is meant to facilitate the execution of the scalpel technique.

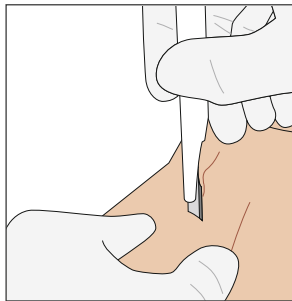


1. Stab

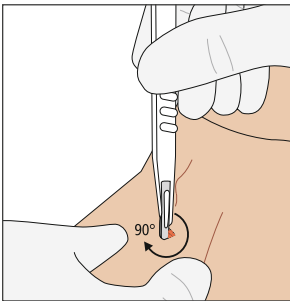

2. Twist

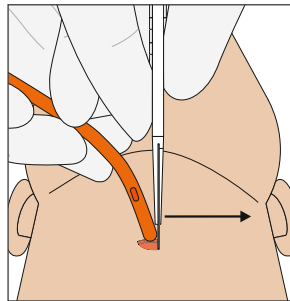


3. Bougie

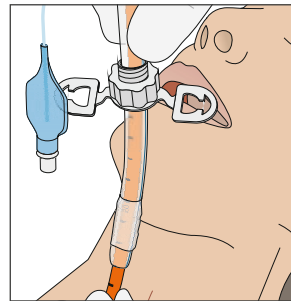


4. Tube

Features

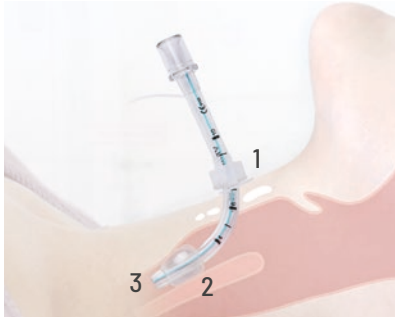

- 1 Adjustable flange
- 2 Thin walled, low pressure cuff
- 3 Smooth tip for ease of insertion

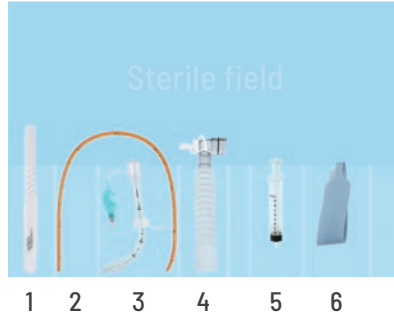

Content

- 1 Scalpel #10
- 2 14 Fr bougie (40 cm)
- 3 6.0 mm I.D. cuffed tube
- 4 Extension tubing
- 5 Syringe 10 ml
- 6 Necktape



The set is packed in a specific sterile pouch, ready to use

Pack size:

24 cm (L) x 13 cm (W) x 4 cm (H)

Unfolded size / sterile area:

56 cm (L) x 39 cm (W)

Order information

ScalpelCric / For single use, sterile

REF	Box
30-08-447-1	1

ScalpelCric / For training, non sterile

REF	Box
30-08-447	1

Additional information



Application Video



Brochure

References

"Anaesthetists must learn a scalpel technique and have regular training to avoid skill fade."

– Baker PA, Weller JM, Greenland KB, Riley RH, Merry AF. Education in airway management. Anaesthesia 2011; 66 (Suppl 2): 101-11

"Training should be repeated at regular intervals to ensure skill retention."

– Difficult Airway Society 2015 guidelines for management of unanticipated difficult intubation in adults.

This product is manufactured without the use of natural rubber latex, unless otherwise specified.

This product does not contain phthalates which require labelling according to CLP Regulation (EC) 1272/2008.

VBM Medizintechnik GmbH

Einsteinstrasse 1 / 72172 Sulz a. N. / Germany

Tel.: +49 7454 9596-0 / Fax: +49 7454 9596-33 / e-mail: sales@vbm-medical.de / www.vbm-medical.de

Connect with us:



CE 0123

FLY004_2.0_EN